
Änderungen vorbehalten
Salvo variazioni
Subject to modification
Tous droits de modification réservés
Reservado el derecho de introducir modificaciones
Wijzigingen voorbehouden

0050GA00F
12.03.2008

ANIMO

Safety and comfort

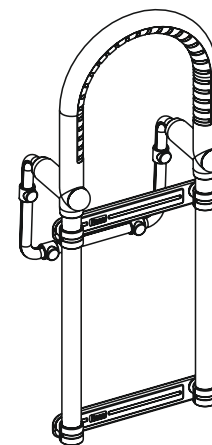PROVEX Industrie GmbH / S.r.l.
Industriezone Nord / 10 Zona Industriale Nord
I-39031 Bruneck / Brunico
Fon: +39 0474 57 15 11
Fax: +39 0474 57 15 71
www.provex.eu
info@provex.eu


PROVEX

made to refresh.

PROVEX

made to refresh.

animo

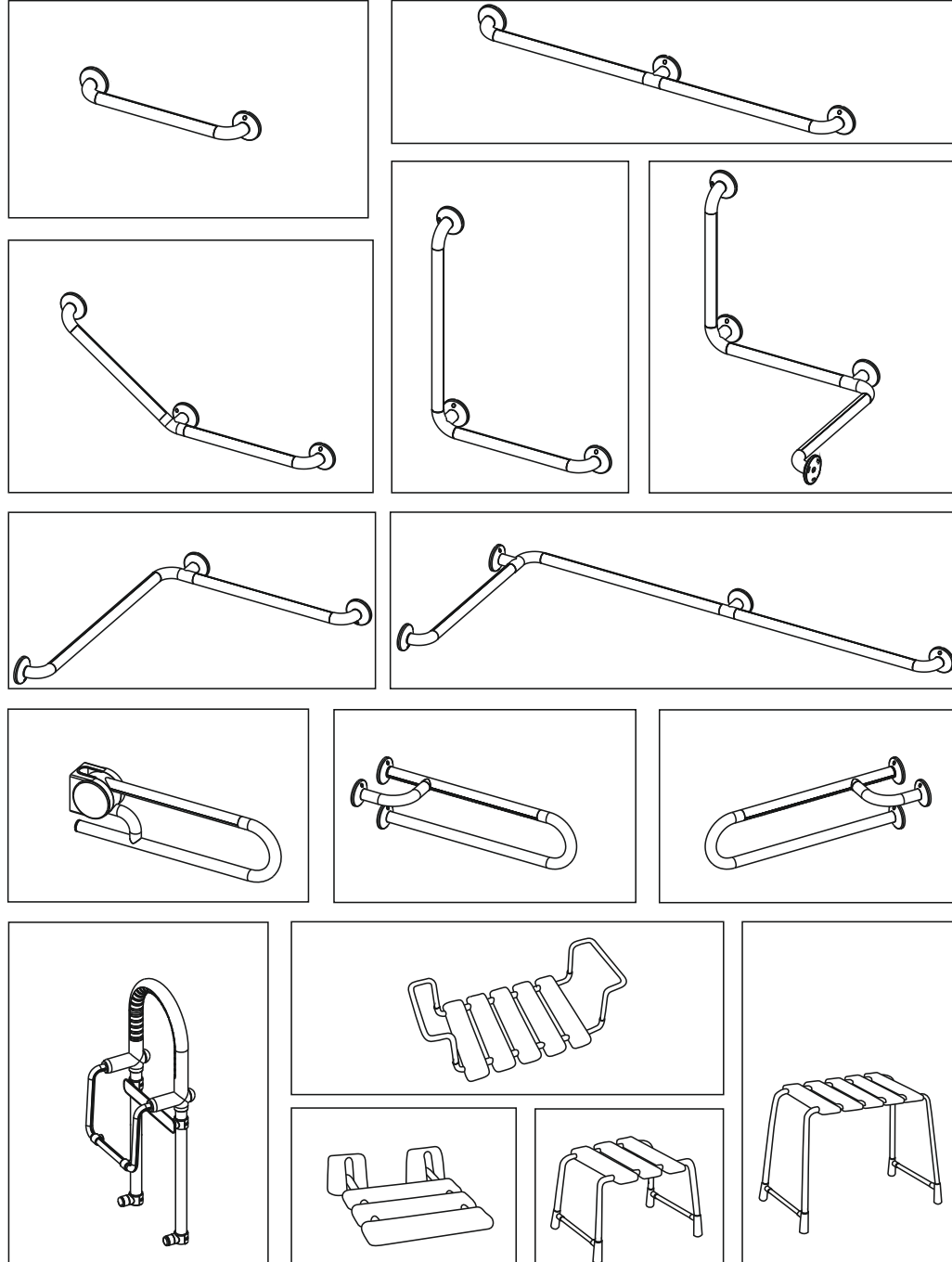Provex - Garantiezusage

Artikelnummer:

Vekäufer:

Kaufdatum:

KA-Nummer:

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Guarantee

The guarantee conditions for this product are as defined by our dealer in the country of sale. Please contact him if necessary. Guarantee and weight capacity only for standard products.

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GB

Cleaning and maintenance

1. Do not use scouring or aggressive detergents, or petrol or thinner.
2. Provex recommends a general Household detergent or "Misterer Propper". After cleaning rinse well with water.
3. Provex recommends for lime deposit a vinegar essence (25%) or household vinegar.
4. Do not use shower gels as there are aggressive to polystyrol syntetic materials.

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GB Mounting with supporting plates (TÜV approved)

- Stick the double adhesive tape (1) on the front side of the plate (2) and pull the protective foil off.
- Fix the plate (2) on the supports (3)
- Insert the rubber stopper (4) on the bracket (5)
- Position the handle on the bath and insert the bracket (5)
- Fix the handle with 2 head cap screws (6), pulling them regularly (do not effort)
- Fix the cover plates (7).

Attention: The Animo safety handle should not be fitted to baths with plastic panels unless adequate support has been provided.

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